



CK SURGITECH

LOCKADO™ CLIP

STRONGER HOLD, MORE PATIENT COMFORT

- Innovative clip design mechanism with defined pre-locking function
- Shortest clip stem in the market
- Additional side teeth on the clips' jaws
- Several different sizes to accommodate clinical needs
- Jaw compression force increased by 15%



NEW



MICRO-TECH™
ENDOSCOPY

YOUR TRUSTED PARTNER IN ENDOSCOPY




LOCKADO™ CLIP

The new LOCKADO™ Clip from MICRO-TECH is the latest development in haemostasis. With the addition of ten teeth that work as anchoring spikes, the user achieves optimal fixation on the tissue of the bleeding area. One of the key features of this clip is its innovative and unique “pre-locking” function. Users can now imitate the clip’s actual closure to confirm correct positioning and haemostasis. Should a change of location be required, the user can open the clip and reposition limitlessly. This option provides many clinical advantages as well as reducing the number of clips used.

MICRO-TECH have also increased the jaw compression force by 15%, enabling a significant improvement in closing pressure. For maximum patient comfort, the clip segment projecting into the lumen was designed to be particularly short. This enables comfortable nutrition and minimizes the risk of injury to the opposing tract. MICRO-TECH offers the LOCKADO™ with different jaw options. Currently available as mini 8mm, standard 11mm and large 16mm, with additional sizes currently under development. LOCKADO™ CLIP brings you the next level of innovation.

SPECIFIC CHARACTERISTICS

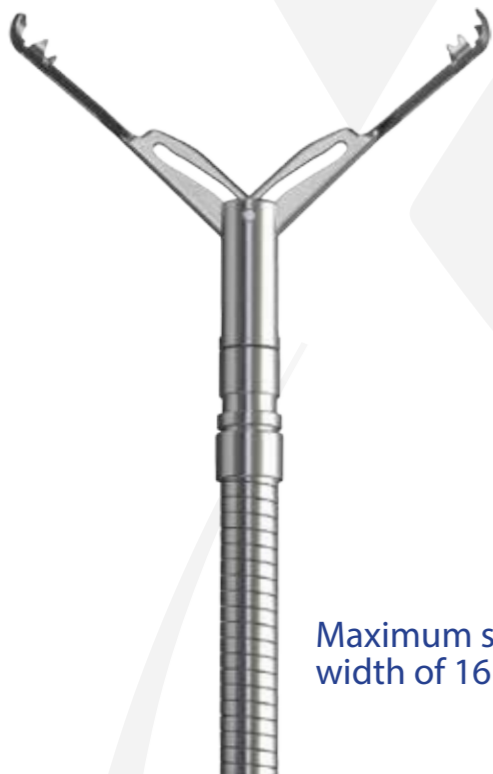

Maximum span width of 16 mm

- Bleeding is better targeted and can be stopped quickly and specifically
- Ten jaw teeth for optimal hold
- Defined pre-locking of the clip arms with maximum contact pressure
- Unlimited repositioning possible prior to fixation
- High fixation due to better adaptation of tissue
- Increase the clips’ compression force by 15%
- Short clip segment for greater patient comfort

SPECIFICATIONS

REF	Clip Max. span width mm	Length mm	Instrument ø mm	Grip colour	Covering
LOCKADO™					
LOCK-C-26-230-C	8	2300	2.6	PINK	PINK
LOCK-D-26-230-C	11	2300	2.6	BLUE	BLUE
LOCK-F-26-230-C	16	2300	2.6	YELLOW	YELLOW

*PACKAGE QTY / 10 PCS

LOCKADO™ Mini

8mm

Lowest jaw profile of 8mm for unobstructed view during placement - shortest stem

LOCKADO™

11mm

For post polypectomy bleeding & defect closure with 11mm jaw - shortest stem

LOCKADO™ Plus

16mm

Widest opening of 16mm for large defect closure with optimal reach



YOUR TRUSTED PARTNER IN ENDOSCOPY

NATIONAL SALES TEAM



Trina Cornwall
Queensland / Territory
M | 0423 036 640



John Lovegrove
QLD / Nth NSW
M | 0412 714 708



Katrina Meyers
NSW & ACT Territory
M | 0406 838 475



Karen Bull
NSW
M | 00401 293 220



Tahlia Atsidaftis
South Australia
M | 0404 045 469



Justine Presnell
Victoria
M | 0452 118 401



Leesa Cutajar
Victoria
M | 0413 204 468



Narelle Fish
WA / TAS Melbourne
M | 0402 656 497

CK Surgitech Pty Ltd
6/14 Technology Drive, Arundel
Queensland Australia 4214

P: 07 5574 5573
E: sales@cksurgitech.com
www.cksurgitech.com

